

Datasheet

HSA-X-S-2G-IN

Ultra High Speed Photoreceiver with InGaAs-PIN Photodiode



The picture shows model HSA-X-S-2G-IN-FST


Features

- **InGaAs-PIN photodiode**
- **Bandwidth 10 kHz – 2 GHz**
- **Amplifier transimpedance gain 5.0×10^3 V/A**
- **Max. conversion gain 4.75×10^3 V/W @ 1550 nm**
- **Spectral range 900 – 1700 nm**
- **Free-space input 1.035"-40 threaded**
- **Fiber optic input available as permanently mounted FC-input**
- **UNC 8-32 and M4 tapped holes for mounting on standard posts with metric and imperial thread**

Applications

- **Spectroscopy**
- **Ultra-fast pulse and transient measurements**
- **Optical triggering**
- **Optical front-end for oscilloscopes and ultra-fast A/D converters**

Block Diagram

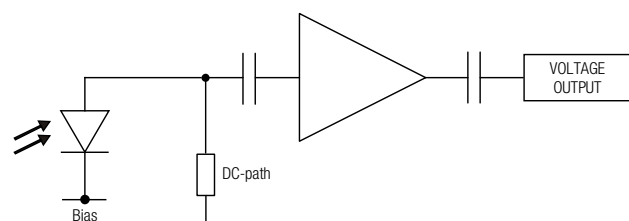

BS01-HSA-X-S_R01

Intended Use

The HSA-X-S-2G-IN photoreceiver consists of an InGaAs photodiode and a subsequent low-noise fixed gain amplifier. It is designed for ultra-fast conversion of small optical signals into equivalent output voltages. Operation is mostly self-explanatory. If in doubt, consult this document or contact support@femto.de.

For safe operation, please refer to the damage thresholds specified in the "Absolute Maximum Ratings", "Temperature Range" and "Power Supply" sections of this document.

The operating environment must be free of smoke, dust, grease, oil, condensing moisture, and other contaminants that could affect the operation or performance.

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E M T O

HSA-X-S-2G-IN_R13/TH, JMa/06MAY2024

Datasheet
HSA-X-S-2G-IN
**Ultra High Speed Photoreceiver
with InGaAs-PIN Photodiode**
Available Versions
HSA-X-S-2G-IN-FST


1.035"-40 threaded flange with internally threaded coupler ring (outer diameter 30 mm) for free space applications. Compatible with many optical standard accessories and for use with various types of fiber connector adapters.

HSA-X-S-2G-IN-FC


Fix/permanent FC fiber connector for high coupling efficiency and excellent conversion gain accuracy.

Related Models
HSA-X-S-1G4-SI-FST

Si-PIN, $\varnothing$ 0.4 mm, 320 – 1000 nm, 1.4 GHz, free space input, 1.035"-40 threaded flange

HSA-X-S-1G4-SI-FC

Si-PIN, integrated ball lens, 320 – 1000 nm, 1.4 GHz, FC fiber connector (fix/permanent)

HSPR-X-I-1G4-SI-FST

Si-PIN, $\varnothing$ 0.4 mm, 320 – 1000 nm, 1.4 GHz, inverting output, free space input, 1.035"-40 threaded flange

HSPR-X-I-1G4-SI-FC

Si-PIN, integrated ball lens, 320 – 1000 nm, 1.4 GHz, inverting output, FC fiber connector (fix/permanent)

HSPR-X-I-2G-IN-FST

InGaAs-PIN, $\varnothing$ 0.1 mm, 900 – 1700 nm, 2 GHz, inverting output, free space input, 1.035"-40 threaded flange

HSPR-X-I-2G-IN-FC

InGaAs-PIN, integrated ball lens, 900 – 1700 nm, 2 GHz, inverting output, FC fiber connector (fix/permanent)

Available Accessories
PS-15-25-L


Power Supply
Input: 100 – 240 VAC
Output: $\pm$ 15 VDC

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY
F E M T O

HSA-X-S-2G-IN_R13/TH, JMa/06MAY2024

Datasheet
HSA-X-S-2G-IN
**Ultra High Speed Photoreceiver
with InGaAs-PIN Photodiode**

Specifications	Test conditions	$V_s = +15\text{ V}$, $T_A = 25\text{ }^\circ\text{C}$, output load impedance $50\ \Omega$, warm-up 20 minutes (min. 10 minutes recommended)
Gain	Transimpedance gain Conversion gain	$5.0 \times 10^3\text{ V/A}$ (@ output load $50\ \Omega$) $4.75 \times 10^3\text{ V/W typ.}$ (@ 1550 nm, output load $50\ \Omega$)
Frequency Response	Lower cut-off frequency (-3 dB) Upper cut-off frequency (-3 dB)	10 kHz 2 GHz ($\pm 15\%$)
Time Response	Rise/fall time (10 % – 90 %)	180 ps ($\pm 15\%$)
Input	Noise equivalent power (NEP) Optical saturation power	16 pW/ $\sqrt{\text{Hz}}$ (@ 1550 nm, 100 MHz) 200 $\mu\text{W AC}$ (for linear amplification, @ 1550 nm) 10 mW CW (to prevent saturation, @ 1550 nm)
Detector	Detector Active area (FST version) Active area (FC version) Spectral range Max. sensitivity	InGaAs-PIN photodiode $\varnothing 100\ \mu\text{m}$ integrated ball lens, suitable for fibers up to 62.5 μm core diameter 900 – 1700 nm 0.95 A/W typ. (@ 1550 nm)
Output	Output voltage range Output reflection S22 Output impedance Output noise	1.9 V_{PP} (@ $50\ \Omega$ output load) for linear operation and low harmonic distortion -7.3 dB (@ $f < 2.5\text{ GHz}$) 50 Ω (terminate with $50\ \Omega$ load) 3.6 mV RMS (24 mV peak-peak) typ. (@ $50\ \Omega$ load, no signal on detector, measurement bandwidth 4 GHz)
Optical Input Connector	Material FST flange Material FST coupler ring Material FC receptacle	1.4305 stainless steel, nickel-plated 1.4305 stainless steel, glass bead blasted nickel silver
Power Supply	Supply voltage Supply current	+15 V 130 mA (depends on operating conditions, recommended power supply capability min. 200 mA)
Case	Weight Material	133 g (0.29 lbs) HSA-X-S-2G-IN-FST incl. coupler ring 110 g (0.24 lbs) HSA-X-S-2G-IN-FC AlMg4.5Mn, nickel-plated
Temperature Range	Storage temperature Operating temperature	-30 $^\circ\text{C}$... +85 $^\circ\text{C}$ 0 $^\circ\text{C}$... +60 $^\circ\text{C}$
Absolute Maximum Ratings	Optical input power (CW) Power supply voltage	12 mW (averaged) 20 V



SOPHISTICATED TOOLS FOR SIGNAL RECOVERY
F E M T O

HSA-X-S-2G-IN_R13/TH, JMa/06MAY2024

Datasheet

HSA-X-S-2G-IN

Ultra High Speed Photoreceiver with InGaAs-PIN Photodiode

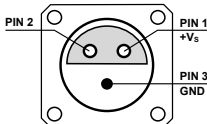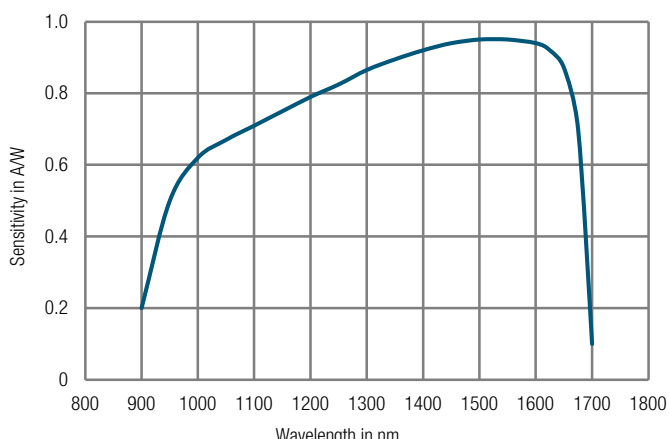Connectors	Input	HSA-X-S-2G-IN-FST	1.035"-40 threaded flange for free space applications and for use with various types of optical standard accessories
		HSA-X-S-2G-IN-FC	FC fiber optic connector (fix/permanent, FC/PC and FC/APC compatible)
	Output	SMA jack (female)	
	Power supply	LEMO® series 1S, 3-pin fixed socket (mating plug type: FFA.1S.303.CLAC52)	
		 Pin 1: +15 V Pin 2: NC Pin 3: GND	
Scope of Delivery	HSA-X-S-2G-IN, internally threaded coupler ring (FST version only), LEMO® 3-pin connector, datasheet, transport package		
Ordering Information	HSA-X-S-2G-IN-FST	1.035"-40 threaded flange for free space applications and for use with various types of optical standard accessories.	
	HSA-X-S-2G-IN-FC	FC fiber optic connector (fix/permanent, FC/PC and FC/APC compatible).	
Spectral Response	 DB-Sens-HSA-X-S-2G-IN_R01		


SOPHISTICATED TOOLS FOR SIGNAL RECOVERY

F E M T O

HSA-X-S-2G-IN_R13/TH, JMa/06MAY2024

Datasheet
HSA-X-S-2G-IN
**Ultra High Speed Photoreceiver
with InGaAs-PIN Photodiode**
Dimensions

HSA-X-S-2G-IN-FST (1.035"-40 threaded free space input)

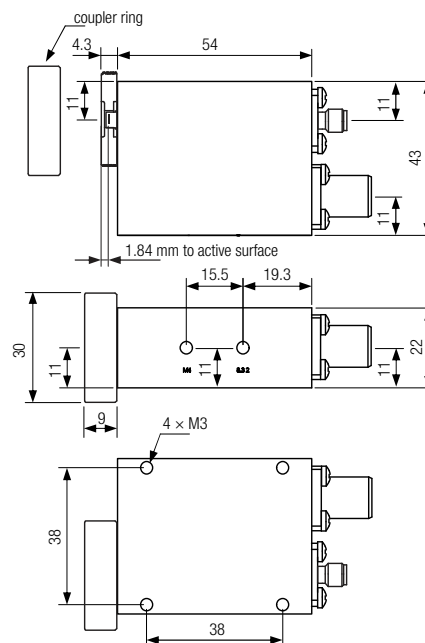

DZ-HS-FST_R1

all dimensions in mm unless otherwise noted


SOPHISTICATED TOOLS FOR SIGNAL RECOVERY
F E M T O

HSA-X-S-2G-IN_R13/TH, JMa/06MAY2024

Datasheet
HSA-X-S-2G-IN
**Ultra High Speed Photoreceiver
with InGaAs-PIN Photodiode**

Dimensions (continued)

HSA-X-S-2G-IN-FC (FC fiber optic connector)

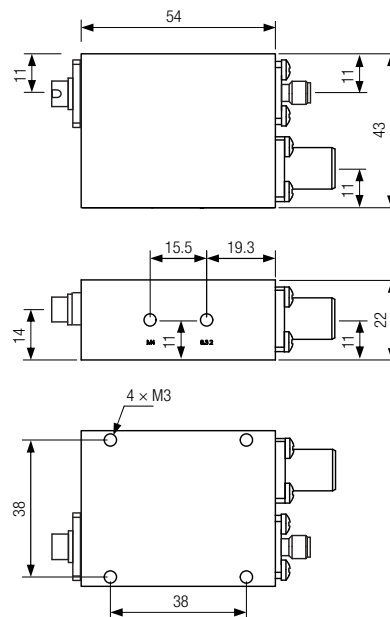

DZ-HS_FC_R1

all dimensions in mm unless otherwise noted

Specifications are subject to change without notice. Information provided herein is believed to be accurate and reliable. However, no responsibility is assumed by FEMTO Messtechnik GmbH for its use, nor for any infringement of patents or other rights of third parties which may result from its use. No license is granted by implication or otherwise under any patent or patent rights of FEMTO Messtechnik GmbH. Product names mentioned may also be trademarks used here for identification purposes only.

© by FEMTO Messtechnik GmbH · Printed in Germany

SOPHISTICATED TOOLS FOR SIGNAL RECOVERY
F E M T O

HSA-X-S-2G-IN_R13/TH, JMa/06MAY2024