


Mirao 52e

Deformable mirror **The looper**

Designed for closed loop Large amplitude Exceptionnal surface quality



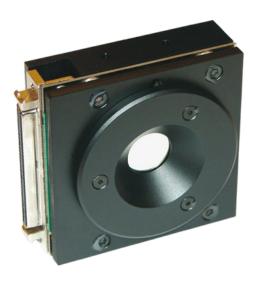

Mirao 52e DATASHEET 2311



Mirao 52 family +


Boost your imaging performance:
Adaptive Optics made easy and efficient.

APPLICATIONS

- + **Ophthalmology :** Explore retinal cells at high resolution (contact our sister company Imagine Eyes for more informations)
- + **Microscopy**: Image deeper in your sample and/or navigate in 3D (for more details visit mu-Imagine website, our division dedicated to microscopy)
- + Quantum applications
- + Beam shaping
- + Laser micromachining
- + Education

FEATURES

- + Fast closed-loop convergence and accurate correction with high linearity and very low hysteresis
- + **Preserved photon budget** with achromatic, highly reflective and continuous membrane
- + Long-term stability with high stabilization option (Mirao 52es), allowing open-loop operation
- + Correction up to 6th Zernike order enabled by 52 electromagnetic actuators
- + **Protected version available** (Mirao 52ep) to prevent membrane mechanical damage

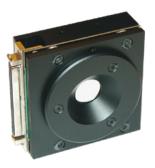

Mirao 52e



Mirao 52ep (protected)



Mirao 52es (stabilized)

Mirao 52e DATASHEET 2311





SPECIFICATIONS*

OPTICAL SPECS

Surface quality Coating Linearity Hysteresis <10 nm RMS Protected silver > 95% < 2%

52 ± 50 μm 15 mm

OPERATING SPECS

Number of actuators Maximum generated wavefront (PV) Effective diameter Spatial frequency correction

Zernike orders up to 6 2.4 ms

Rise time Temporal stability

OPERATING SYSTEM

< 10 nm RMS over 12h (stabilized option)

MISC

Dimension / Weight (Mirao 52e unit only) Dimensions / Weight (Mirao 52e controller) Working environment 64 x 64 x 23 mm³ / 490 g 24 x 23 x 10 cm³ / 3 kg 20-25°C, 20-80% RH USB 2.0 / 50 W


Interface / Power consumption


Windows 10

*Subject to changes without further notice

Zernike Order (PV) wavefront generation 1 2 ±50μm ±50μm ±50μm 3 ±25μm ±10μm ±10μm ±25μm 4 ±15μm ±8μm ±8μm ±8μm ±15μm

DIMENSIONS** (mm)

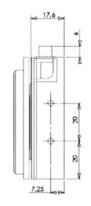

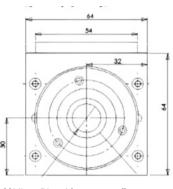

**Mirao 52e without controller

Mirao 52e DATASHEET 2311



SOFTWARE

WAVETUNE

WAVETUNE is a unique software that seamlessly combines wavefront measurement and correction features with extensive instrument diagnostics.

This software contains all the necessary tools to calibrate the Deformable Mirror (DM). It can also operate the DM in closed-loop with HASO wavefront sensor, as well as in open-loop and perform beam shaping.



WAVEKIT BIO

WAVEKIT BIO is a Software Development Kit (SDK), available in C++ and Python, specifically designed for microscopy applications. In particular, it contains all the necessary functions to implement sensorless AO, using image-based iterative algorithms (e.g. 3N).



MOUNTING & ACCESSORIES

Several mounting options are available, including adaptors for the most common mechanical stages, to simplify integration of any Mirao 52 device into an optical setup.

